

Technical documentation

Product name:	qRE Filipendula ulmaria (L.) Maxim., flowering tops
Substance:	Filipendula ulmaria (L.) Maxim., flowering tops dry extract
Plant source common names:	en: Meadowsweet ; fr: Reine des prés
Reference:	E0067
Packaging:	100 mg in a 1.5 ml borosilicate amber vial
Storage conditions:	Keep container closed. Protect from light and moisture. Keep at -18 °C
Retest:	12 months

Botanical identification of plant source

Plants in our botanical garden are identified and a herbal voucher is prepared by an expert botanist.

Each batch collected for extraction is verified and identified.

Reference: Flora Europaea, Cambridge University Press, 1968, Vol 2, p 6

Method of production of dry extract

Whole plant or plant parts are collected, freeze-dried and coarsely ground. Extraction is performed by maceration in 50 % (v/v) aqueous ethanol for 48 hours at room temperature. Ethanol is then evaporated under reduced pressure at less than 40 °C and the aqueous residue is freeze-dried.

Organoleptic characteristics of dry extract

Colour: Yellow beige

Odour: Characteristic

Form: Fine powder

Recommended methods for use

Weight a precise weight of qRE and solubilise in the recommended solvent at the concentration indicated in the HPLC or HPTLC method described in this document.

Sonicate for 90 seconds (70 W). Filter on a 0.45 µm PVDF membrane and put the resulting solution into HPLC dispenser or apply on the HPTLC plate.

Dose and analyse your extract with qRExtract using the HPLC / HPTLC methods described in this document or using your own methods.

HPTLC

Detection of hyperoside, quercetin, kaempferol and rutin

Layer: 10 × 10 cm HPTLC Nano-Sil-20 UV 254 (Carl Roth ref. N084.1)

Thin layer conditionnement: 1 h at room temperature and 33 % relative humidity

Elution solvent:	<u>Elution solvent compound</u>	<u>Volume (ml)</u>
	ethyl acetate	15
	formic acid	1
	H ₂ O	1

Developing distance: 70 mm from the lower edge

Initial spot volume and concentration:

hyperoside:	2 µl of a 0.02 % (w/v) solution in methanol
quercetin:	2 µl of a 0.02 % (w/v) solution in methanol
qRE:	5 µl of a 2.4 % (w/v) solution in 50 % (v/v) aqueous ethanol
kaempferol:	2 µl of a 0.02 % (w/v) solution in 50 % (v/v) aqueous ethanol
rutin:	1 µl of a 0.1 % (w/v) solution in methanol

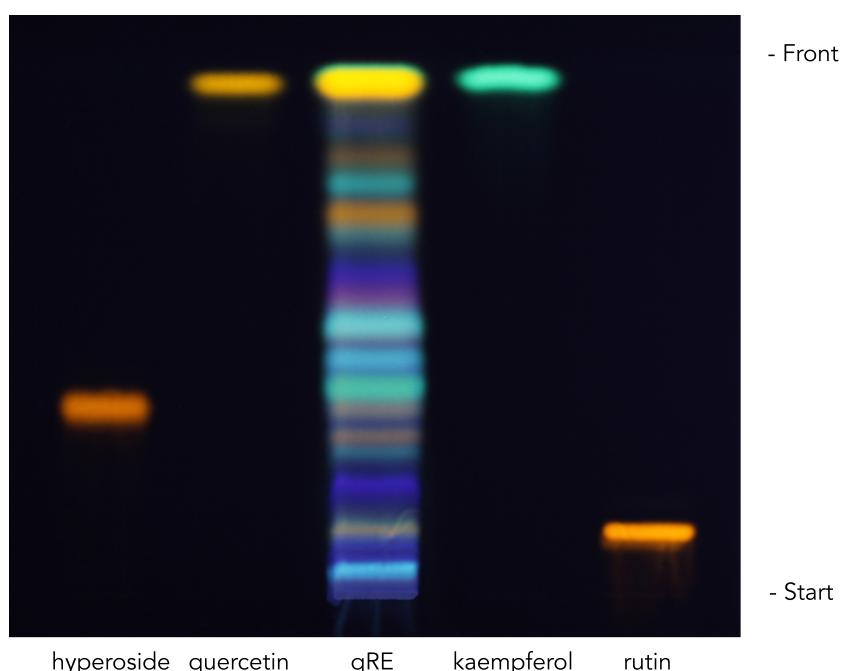
Reagent mixture:

Natural products - polyethylene glycol reagent (NP/PEG)

Preparation: Dissolve 0.25 g of diphenylboric acid 2-aminoethylester and 1.25 g of polyethylene glycol 400 in 25 mL of methanol.

Dip the plate in the reagent mixture and dry for 15 minutes at room temperature.

Expose to UV light at 365 nm.



HPTLC

2nd detection of quercetin and kaempferol

Layer: 10 × 10 cm HPTLC Nano-Sil-20 UV 254 (Carl Roth ref. N084.1)

Thin layer conditionnement: 1 h at room temperature and 33 % relative humidity

Elution solvent:	<u>Elution solvent compound</u>	<u>Volume (ml)</u>
	toluene	10
	ethyl acetate	8
	formic acid	2


Developing distance: 70 mm from the lower edge

Initial spot volume and concentration:

- | | |
|-------------|---------------------------------------------------------------|
| quercetin: | 2 µl of a 0.02 % (w/v) solution in methanol |
| qRE: | 4 µl of a 2.4 % (w/v) solution in 50 % (v/v) aqueous ethanol |
| kaempferol: | 2 µl of a 0.02 % (w/v) solution in 50 % (v/v) aqueous ethanol |

Reagent mixture: Natural products - polyethylene glycol reagent (NP/PEG)

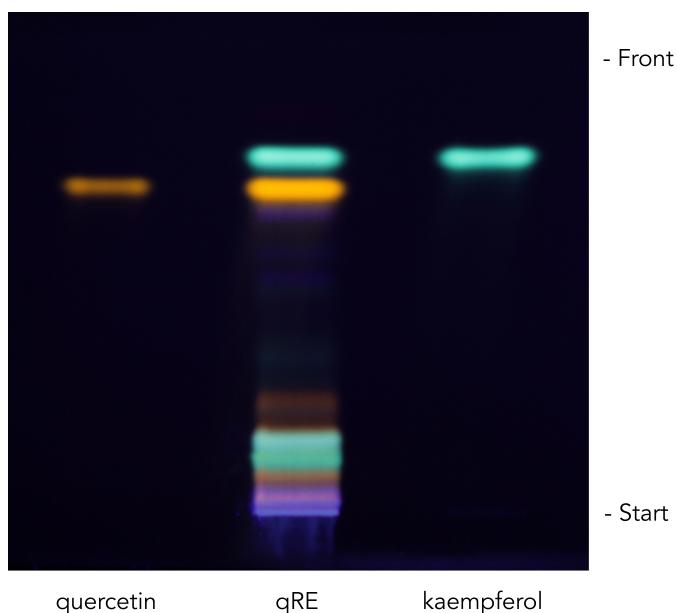Preparation: Dissolve 0.25 g of diphenylboric acid 2-aminoethylester and 1.25 g of polyethylene glycol 400 in 25 mL of methanol.
Dip the plate in the reagent mixture and dry for 15 minutes at room temperature.
Expose to UV light at 365 nm.



HPLC

Precolumn: Ascentis® Express C18 0.5 cm × 3.0 mm 2.7 µm

Column: Ascentis® Express C18 15 cm × 3.0 mm 2.7 µm

Sample: 8 µl 0.5 % qRE (w/v) solution in 25 % (v/v) aqueous ethanol

Flow: 0.45 ml/min

Temperature: 25 °C




Mobile phase: A: 0.1 % formic acid (v/v) in water
B: 0.1 % formic acid (v/v) in acetonitrile

Detection: Diode Array Detector, 280 nm

Gradient:

Time (mn)	A %	B %
0	97	3
30	90	10
100	83	17
110	76	24
130	0	100
140	0	100

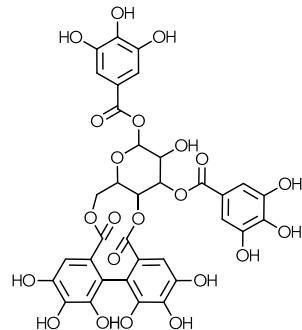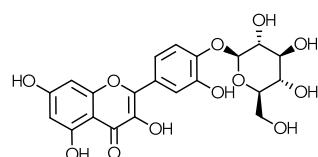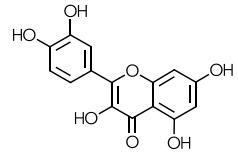Quantified substances

Compound	CAS No	2D Structure	Peak No
Digalloyl-HHDP-glucose isomer	NA		1, 2
Quercetin 4'-O-glucoside	NA		8
Unknown	NA	NA	3, 4, 5, 6, 7, 9, 10
Quercetin	117-39-5		11
Kaempferol	520-18-3		12